
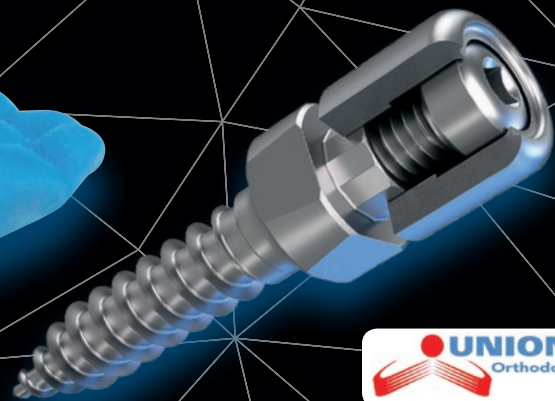
BENEFIT[®] .KFO

EASY DRIVER[®]
THE GUIDING SYSTEM FOR TADS

Digital skeletal anchorage
Digitale skeletale Verankerung

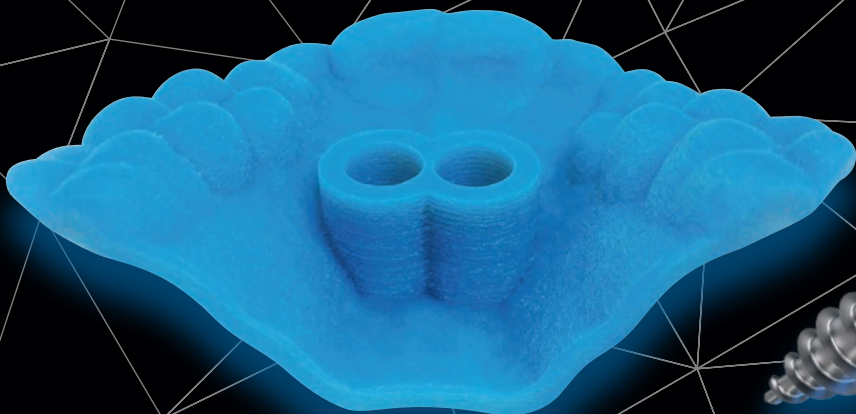www.psm.ms – implants for professionals

Digital Technology

– EASY DRIVER

Digitale Technologie

– EASY DRIVER

A project developed in cooperation with Uniontech Laboratories (Parma/Milano), the FACE Group and Dr. Renzo di Gabriele.

The worldwide leading **BENEFIT®** Miniscrew System offers numerous advantages for skeletal anchorage. It allows an easy insertion of the miniscrew by the orthodontist. Additionally it overcomes limits of traditional appliances and allows tooth movement previously impossible to achieve.


Difficult anatomical situations in the palate can be challenging for evaluation of the ideal implant position and for the fabrication of the appliance.

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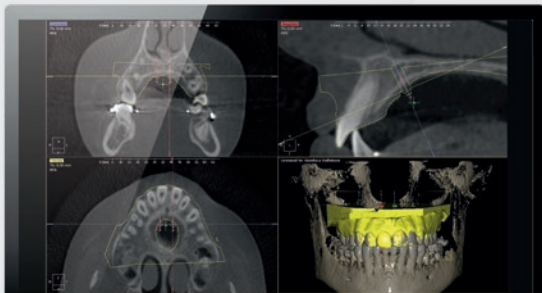
Transfer of the "CBCT" data to evaluate the ideal implant position with virtual **BENEFIT®** Screws.

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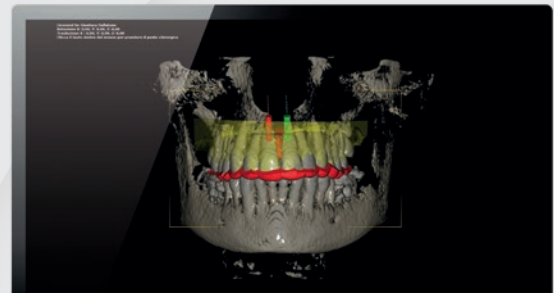
Visualization of the gauge in place.

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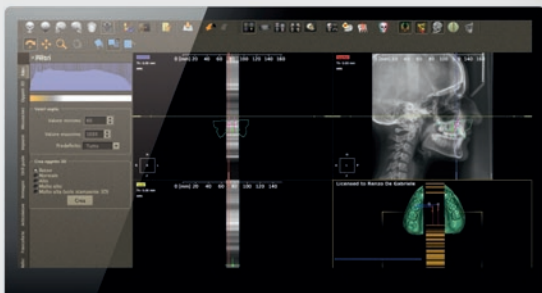
The latero-lateral x-rays in combination with the digital software result in a precise display of the anatomy of the palate.

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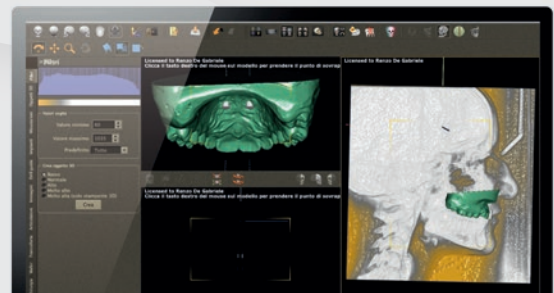Three dimensional verification showing position, inclination and distance of the implants in place.

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The newly developed technique called **„3DRM – resource matching“** allows

- Perfect positioning of the **BENEFIT®** miniscrews for fabrication of the **EASY DRIVER** gauge
- Fabrication the gauge from a first print

Protocoll **BENEFIT®** – **EASY DRIVER**

- Impression taking
- Digitalization of the model
- Adaptation of the model with CBCT or latero-lateral x-ray
- Positioning of the **BENEFIT®** miniscrews
- Re-confirmation or correction of data by the Orthodontist
- Fabrication of gauge and appliance

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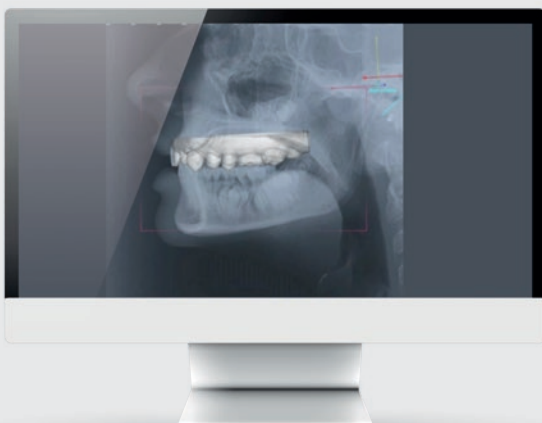
benefit

EASY DRIVER

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- Considerably higher success rate
- Best possible stability
- Secure and easy placement with drill and screw guides
- Planning of ideal screw size and position
- Reduced planning and treatment time
- Implants and appliance placement in one visit

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A new innovative method for the secure insertion of the **BENEFIT®** miniscrews

The surgical **EASY DRIVER** gauge is manufactured with the AM technique (Additive Manufacturing). It allows the Orthodontist to place **BENEFIT®** miniscrews with precise guidance to achieve the perfect angulation and parallelism.

Implantation and appliance placement in one session


The laboratory produces the gauge and also the requested appliance from the digital files. The Orthodontist receives gauge and finished appliance together for final implantation.

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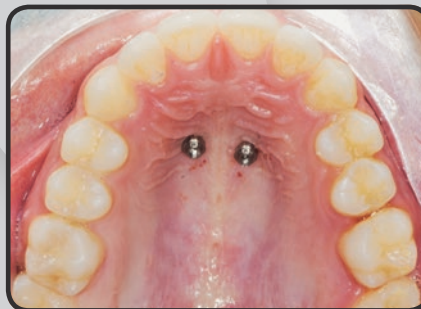Implantation Planning Software Software zur Implantations- Planung

The latest Treatment Planning Software "**EASY DRIVER**" can be obtained under e-mail info@uniontech.it.

The cost of the software is 300 EUR (net, subject to VAT if applicable) including all future updates. No annual fees or additional costs.

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
BENEslider with lingual brackets.
BENEslider mit Lingualbrackets.

Some possible Indications

- Maxillary molar distalization and mesialization
- Molar uprighting
- Anchorage for upper anterior teeth and molars
- Temporary tooth replacement
- Palatal expansion
- Indirect anchorage

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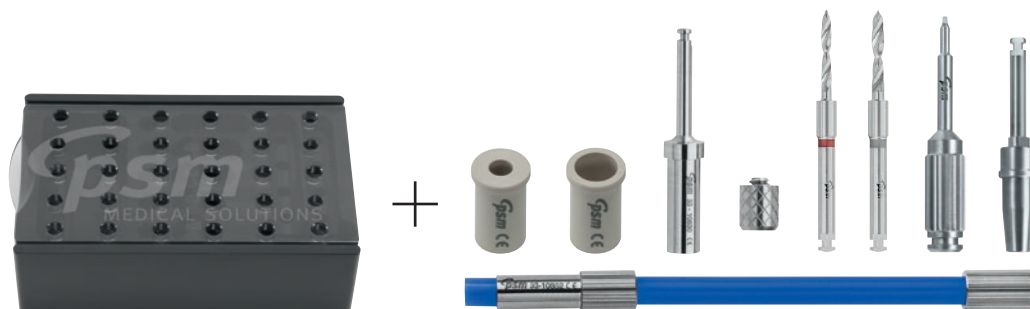**EASY DRIVER
COMPLETE SET
EASY DRIVER
KOMPLETT-SET**

Recommendation
Empfehlung

BENEsider Video

33-56000	EASY DRIVER Complete Set consisting of:	EASY DRIVER Komplett-Set bestehend aus:
	Storage	Aufbewahrung
73-31995	EASY DRIVER Tray black, empty	EASY DRIVER Behälter schwarz, leer
	Instruments	Instrumente
33-10834	Drill guide sleeves Peek EASY DRIVER, 2 ea.	Bohrhülsen Peek EASY DRIVER, 2 Stk.
33-10835	Screw guide sleeves Peek EASY DRIVER, 2 ea.	Schraubhülsen Peek EASY DRIVER, 2 Stk.
33-10830	BENEFIT Screwholder EASY DRIVER	BENEFIT Schraubenhalter EASY DRIVER
33-10850	Guiding cap EASY DRIVER, 2 ea.	Führungskappe EASY DRIVER, 2 Stk.
33-10852	Cap Remover EASY DRIVER	Kappenentferner EASY DRIVER
33-10814	EASY DRIVER Drill 1.4 mm, Dental, red	EASY DRIVER Vorbohrer 1,4 mm, Dental, rot
33-10818	EASY DRIVER Drill 1.8 mm, Dental, grey	EASY DRIVER Vorbohrer 1,8 mm, Dental, grau
33-18288	Laboratory Finger Screwdriver	Labor Finger Schraubendreher
33-18244	Extension, +16 mm Dental	Verlängerung, +16 mm Dental

**EASY DRIVER
BASIC SET
EASY DRIVER
BASIS-SET**

33-56005	EASY DRIVER Basic Set consisting of:	EASY DRIVER Basis-Set bestehend aus:
	Storage	Aufbewahrung
73-31995	EASY DRIVER Tray black, empty	EASY DRIVER Behälter schwarz, leer
	Instruments	Instrumente
33-10834	Drill guide sleeves Peek EASY DRIVER, 2 ea.	Bohrhülsen Peek EASY DRIVER, 2 Stk.
33-10835	Screw guide sleeves Peek EASY DRIVER, 2 ea.	Schraubhülsen Peek EASY DRIVER, 2 Stk.
33-10830	BENEFIT Screwholder EASY DRIVER	BENEFIT Schraubenhalter EASY DRIVER
33-10850	Guiding cap EASY DRIVER, 2 ea.	Führungskappe EASY DRIVER, 2 Stk.
33-10852	Cap Remover EASY DRIVER	Kappenentferner EASY DRIVER
33-10814	EASY DRIVER Drill 1.4 mm, Dental, red	EASY DRIVER Vorbohrer 1,4 mm, Dental, rot

for 2.0 mm screws
für 2,0 mm Schrauben

33-10814 DENTAL

EASY DRIVER Drill 1.4 mm, Dental, red
EASY DRIVER Vorbohrer 1,4 mm, Dental, rot

for 2.3 mm screws
für 2,3 mm Schrauben

33-10818 DENTAL

EASY DRIVER Drill 1.8 mm, Dental, grey
EASY DRIVER Vorbohrer 1,8 mm, Dental, grau

**DRILLS
BOHRER**

All drills:
max. 800 rpm
Alle Bohrer:
max. 800 U/min

Packaging: 1ea.
VPE: 1Stk.

33-10834

Drill guide sleeves Peek EASY DRIVER, 2 ea.
Bohrhülsen Peek EASY DRIVER, 2 Stk.

33-10835

Screw guide sleeves Peek EASY DRIVER, 2 ea.
Schraubhülsen Peek EASY DRIVER, 2 Stk.

33-10852

Cap Remover EASY DRIVER
Kappenentferner EASY DRIVER

33-10830

BENEFIT® Screwholder for EASY DRIVER System
BENEFIT® Schraubenhalter für das EASY DRIVER System

33-10850

BENEFIT® Guiding cap EASY DRIVER, 2 ea.
BENEFIT® Führungskappe EASY DRIVER, 2 Stk.

optional

33-18288


Laboratory Finger Screwdriver, fixed
Labor Finger Schraubendreher, fest

optional

33-18244 **QB** DENTAL

Extension, Dental, +16 mm
Verlängerung, Dental, + 16 mm

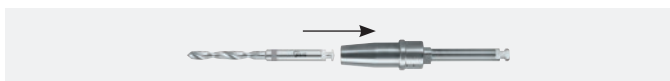The extension of 16 mm works for all rotary instruments of the **quattro®** and **BENEFIT®**-Systems with a latch type DENTAL attachment (max. rpm 600).

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Phone +49 (4402) 8 63 78-0 | Fax +49 (4402) 8 63 78-99

E-Mail info@RKSortho.de

Specifications and design are subject to change without prior notice. All pictures, photos and product descriptions are the property of **psm** MEDICAL SOLUTIONS and UNIONTECH ORTHODONTIC LAB. Utilization and copies by third parties have to be authorized by **psm** MEDICAL SOLUTIONS. All rights reserved!

All trademarks and brands mentioned in this brochure are subject to unrestricted regulations of the trade mark laws valid and right of ownership of the registered owner.

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RKSortho GmbH